

M5 GFP-Tag(7G9) mAb 使用说明书

产品名称	单位	货号
M5 GFP-Tag (7G9) mAb	50 μ l	MF081-T
M5 GFP-Tag (7G9) mAb	100 μ l	MF081-01

【STORAGE】

Store at -20°C. Stable for one year from the date of shipment.

【BACKGROUND】

The green fluorescent protein (GFP) was originally identified as a protein involved in the bioluminescence of the jellyfish *Aequorea victoria*. GFP cDNA produces a fluorescent product when expressed in prokaryotic cells, without the need for exogenous substrates or cofactors, making GFP a useful tool for monitoring gene expression and protein localization in vivo. Several GFP mutants have been developed, including EGFP, which fluoresce more intensely than the wildtype GFP and have shifted excitation maxima, making them useful for FACS and fluorescence microscopy as well as double-labeling applications. GFP is widely used in expression vectors as a fusion protein tag, allowing expression and monitoring of heterologous proteins fused to GFP.

【REFERENCES】

- 1.Prasher, D.C., et al. 1992. Primary structure of the *Aequorea victoria* green fluorescent protein. *Gene* 111: 229-233.
- 2.Chalfie, M., et al. 1994. Green fluorescent protein as a marker for gene expression. *Science* 263: 802-805.

【SOURCE】

This monoclonal antibody is produced by immunizing animals with full length recombinant GFP.

【SPECIFICITY】

GFP-Tag (7G9) Mouse mAb detects GFP and GFP fusion proteins.

【REACTIVITY】

All

【ISOTYPE】

Mouse IgG1

IMPORTANT


Use an anti-MOUSE secondary antibody to detect the 7G9 antibody.

RECOMMENDED ANTIBODY DILUTIONS

Western blotting	1:5000
Immunofluorescence	1:500
Immunoprecipitation	1:50
ELISA	1:2000


For Western blots, incubate membrane with diluted antibody in 5% w/v nonfat dry milk, 1X tbs, 0.05% Tween-20 at 4°C with gentle shaking, overnight.

APPLICATION DATA

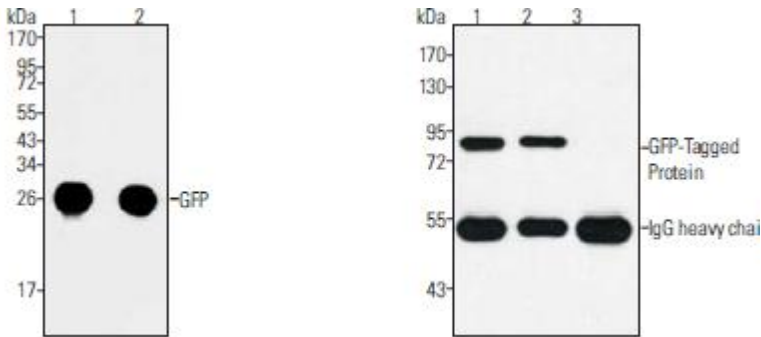

Western blot analysis of over-expressed GFP-tagged protein in 293T cell lysate, using M5 GFP-tag (7G9) Mouse mAb. The antibody dilutions are 1:5000 (lane 1) and 1:10000 (lane 2). Each lane was loaded with 10µg of cell lysate.

IP of extracts from 293T cells transfected (lane 1, 2) with a GFP-tagged protein and a mock transfection (lane 3, using the same protein without the GFP tag), using M5 GFP-Tag (7G9) mouse mAb and probed on Western blot using the same antibody. Dilution: 1:50 (1), 1:100 (2), 1:50 (3).

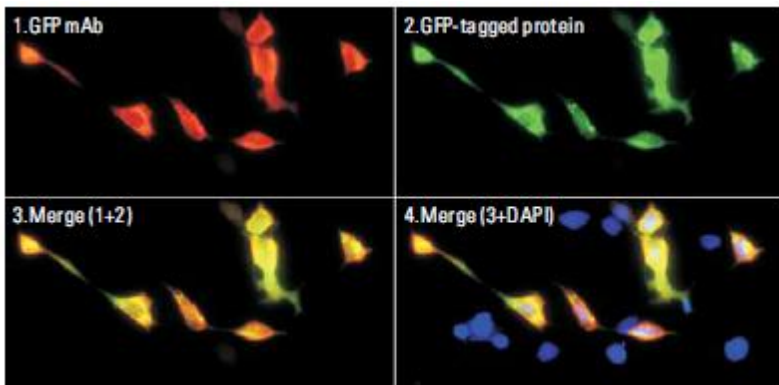

IF analysis of 293T cells transfected with a GFP-tagged protein, using M5 GFP-Tag (7G9) Mouse mAb at a 1:500 dilution.

【备注】

本产品仅供科研使用。在确认产品质量出现问题时，本公司承诺为客户免费更换等量的质量合格产品。